

MAGIC DISC UBS 315


Modulated and synchronized photon, magnetic and scalar therapy in the most compact form. This therapy field applied simultaneously is profoundly more efficient as when it is applied sequential; all living beings will benefit from such a therapy. This effect is remarkable in detoxification, improvement of body fluid rheology, metabolic stimulation and correction, cellular and matrix regeneration.

MAIN FEATURES:





- UBS 315 comes with an intensive central directed photon source; this source emits white light at a temperature of 8000 Kelvin and, alternative, with red light at 633 nm at 6000 mc; both lights are known for many positive health benefits such as: stimulates pineal gland and thus the melatonin secretion for a restful sleep. Red light stimulates collagen and restoration of connective tissue and skin, speeds up wound healing
- Photonic beam is always modulated with the solar noise identical to the sun light spectrum; “white sun light” is the most important frequency to create structured water clusters. Structured water enters into resonance with biological cells and triggers self-healing forces
- Synchronically with the light, the disc emits a pulsating magnetic field, created by a Tesla coil; this field comes with minimum 40% scalar waves, the information therapy. The magnetic field is modulated by Schuman frequency 7.8 -8Hz. Tesla coil field is 150 micro Tesla
- It is possible to modulate additionally through external connected audio signals (do not use mobile phone or tablets); in case of an external audio modulation, both the photons beam and the magnetic field are modulated with the external signals
- Remedies, water, water soluble fluids, homeopathic carriers can be imprinted, or bio-signals can be transmitted directly to the body
- UBS 315 comes with a miniature extra powerful and efficient rechargeable lithium battery 3.6V and 350 mAh lasting for a continuous usage of 2 hours; rechargeable through USB cable
- Main unit sizes: Día 95 mm x 11 m height, Weight 60 grams
- If you have a NSA-7 you can beautifully modulate your UBS 315 light with additional health benefit

SCOPE OF DELIVERY

- Magic therapy disc UBS 315 control unit
- Mini USB cable (charging)
- Audio signal connecting cable
- Mobile holder (ULD Mobile)
- Tesla Transmitter (UBS 315 Trans)

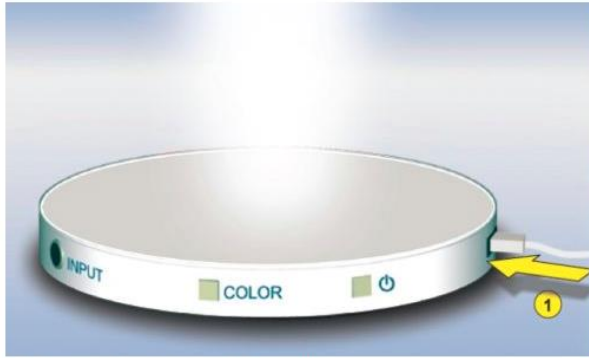

UBS 315 Photon, magnetic scalar synchronic therapy

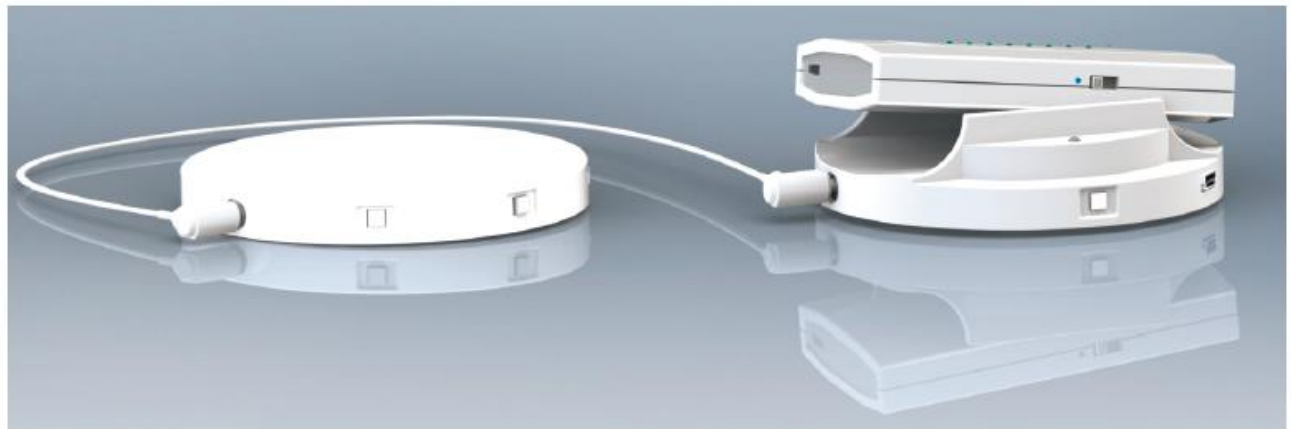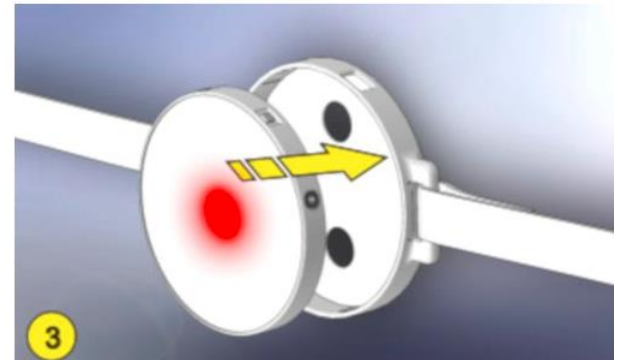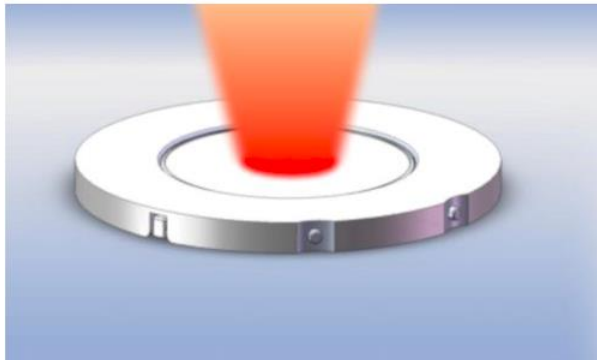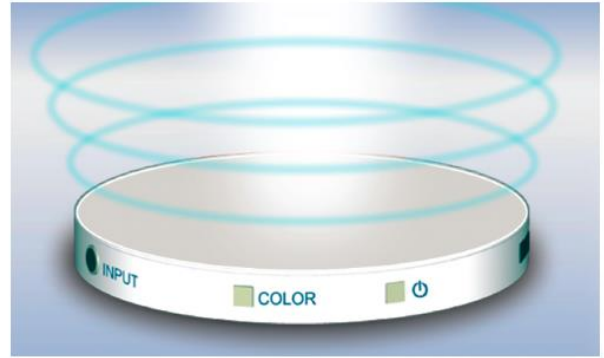


Brain waves modulation with NSA-07

HOW IT WORKS:

The LED diodes are stimulating the water molecules inside the body (over 76%), structuring the water in hexagonal clusters (jelly water), which go into resonance with the cells. Sick cells are thus stimulated to detoxify, relieve pain and regress inflammation, boost cellular communication, repair and restore. This therapy can be applied directly on the affected body zones as well as by imprinting water and remedies which are consumed thereafter.

The magnetic field generated by a Tesla coil contains over 40% scalar waves, the frequencies carrying information; this information interacts with the DNA and microtubules, correcting and restructuring the sick body parts with healing and beneficial information. The synergy

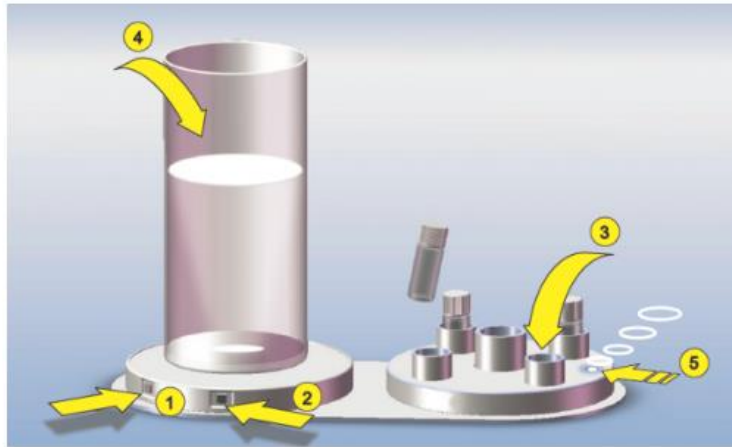

UBS 315 IS SIMPLE AND ROBUST, EASY TO OPERATE AND TO APPLY THE THERAPY

Uncountable health benefits are a result either of a direct therapy applied to the body or of drinking the imprinted water or of orally taking the customized generated remedies; the health benefits will become remarkable very quickly:

- Cellular regeneration and deep detoxification
- Genetic modulation and genetic correction
- Pain relief, restoration of connective tissue, cartilages and energizing of muscular skeletal system
- Supportive therapy of all multiple chronic conditions
- Improvement of rheological properties of body fluids, blood, lymph
- Mitochondrial activation, energy boost
- Stress relief, depression, anxiety and insomnia therapy, emotional trauma cure
- Hormonal, metabolic and self-healing stimulation
- Cellular regeneration especially after injuries, surgeries, cell loss
- Anti-inflammatory
- Energizing, stamina, boost performance, clear mind
- Brain stimulation and hemispheres synchronization

HOW TO APPLY THE THERAPY:

- ***Nervous system:*** Solar plexus, Back of the head, Vertebra C7
- ***Immune system:*** Thymus, On and below the navel, Back of the head
- ***Endocrine system:*** Top of the skull, thyroid, pubic region, center of the spine
- ***Respiratory system:*** on the middle of sternum, under the collar bones, left and right joints shoulder thorax, on the back side between shoulder blades
- ***Circulatory system:*** on each kidney, underneath vertebra C7, back of the head

- **Detox and Lymphatic system:** om vertebra C7, under collar bones, back of the head, liver, each kidney
- **Gastro-intestinal tract:** Solar plexus, un the navel and below navel 1 inch, left and right of the knees, back of the head, underneath should blades left and right
- **Kidneys and Urinary Tract:** back of the head, each kidney, inside the legs on the ankles, bladder and along the spine left and right (bladder meridian)
- **Liver and Gallbladder:** underneath right ribcage, back of the head, right shoulder
- **Joints and spine:** treat always the left and the right knees, shoulders, elbows, feet, hand wrists; treat inside and outside of the joints, on the top and underneath the joints; for the spine insist on vertebra C7 and sacrum, go along the spine up and down, always left and right of the spine not on the spine directly
- **Jaws, teeth, mouth and gum:** on the sinuses forehead and left and right the nose, underneath jaws, back head, behind ears and on temporal lobes
- **Eyes:** use only the red LED light, never use the white radiation: close the eyes and apply on each eye, underneath eyes, back of the head, temporal lobes, forehead
- **Pain:** apply directly on the painful zone(s) or on pain trigger points indicated by your health care fiver
- **Wounds/Injuries:** Never apply directly on the wounds, burns, only around them
- **Skin eczemas:** Never apply directly on the wounds, burns, only around them
- **You can apply UBS 315 on acupuncture points, reflex zones**

MAINTENANCE AND SAFETY:

1. Both LED diodes for white and red light have a lifespan of 20000 hours, corresponding to a usage of 4 hours per day, it means the device is lasting 10 years
2. The device enclosure is biological compatible; never place the device in humid spaces or under the direct sunlight, do not put objects heavier than 0.5 grams
3. Avoid that the small windows where light is being radiated are getting dirty; clean the enclosure with a soft cotton cloth slightly damp, do not use solvents; pay attention that no water is penetrating inside of the enclosure
4. Always charge the battery, never store the device with empty battery; charging time is 2-4 hours